


Participant Instructions for Viewing Precipoint Scanned Images for The Leukaemia Diagnostic Interpretation (LDI) Programme

PLEASE NOTE – The images for the LDI programme are now being hosted by PreciPoint and will therefore appear differently to previous exercises. The image will also be embedded on the data entry page so that participants can view it directly from there without being directed to PreciPoint's website.


Once on the image webpage, the image is accessible using the controls on the image pane.




1. There are several ways to zoom in on the image.
 - a. The times magnification at the top left of the screen

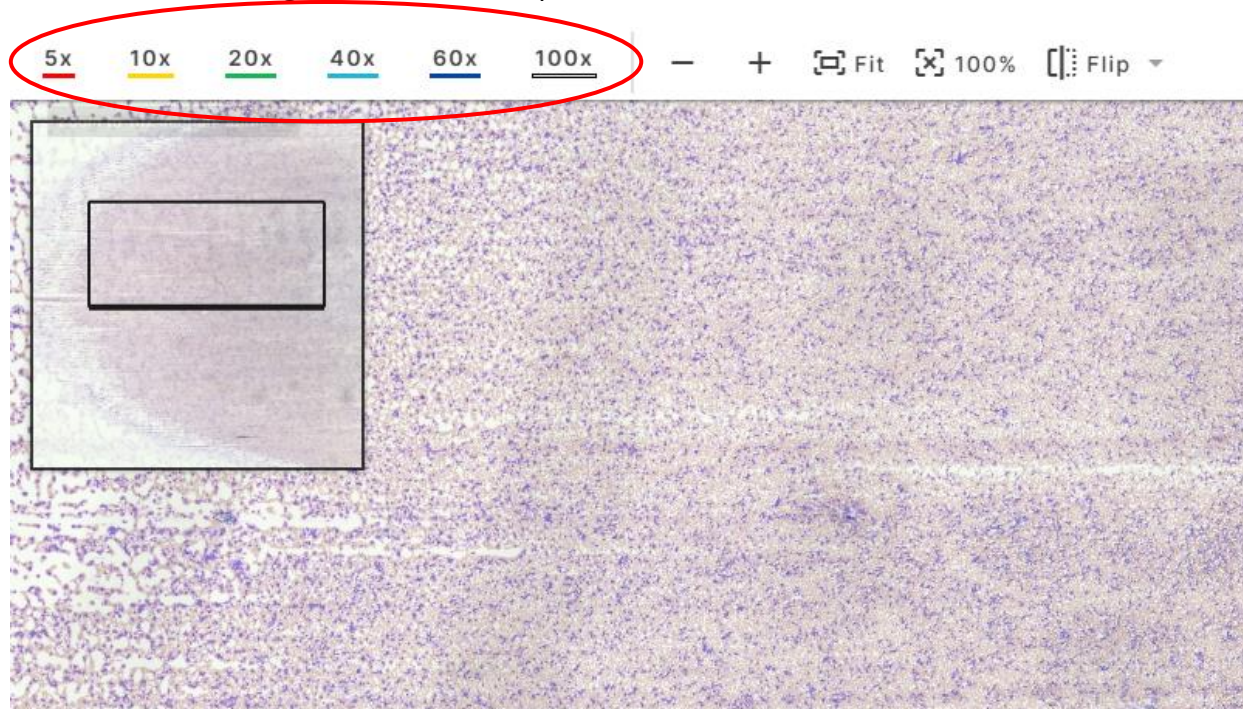

b. The – and +

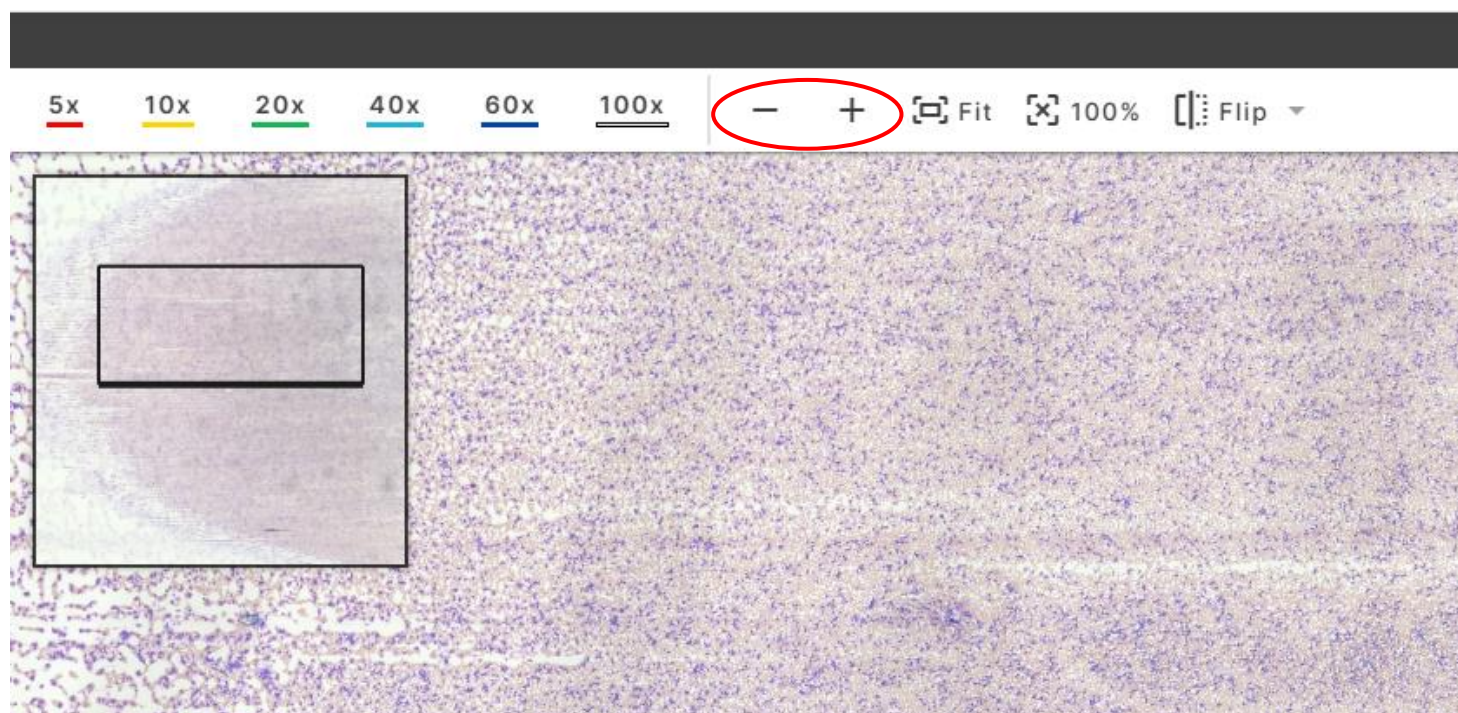

c. The FIT icon will return the magnification to fully zoomed out

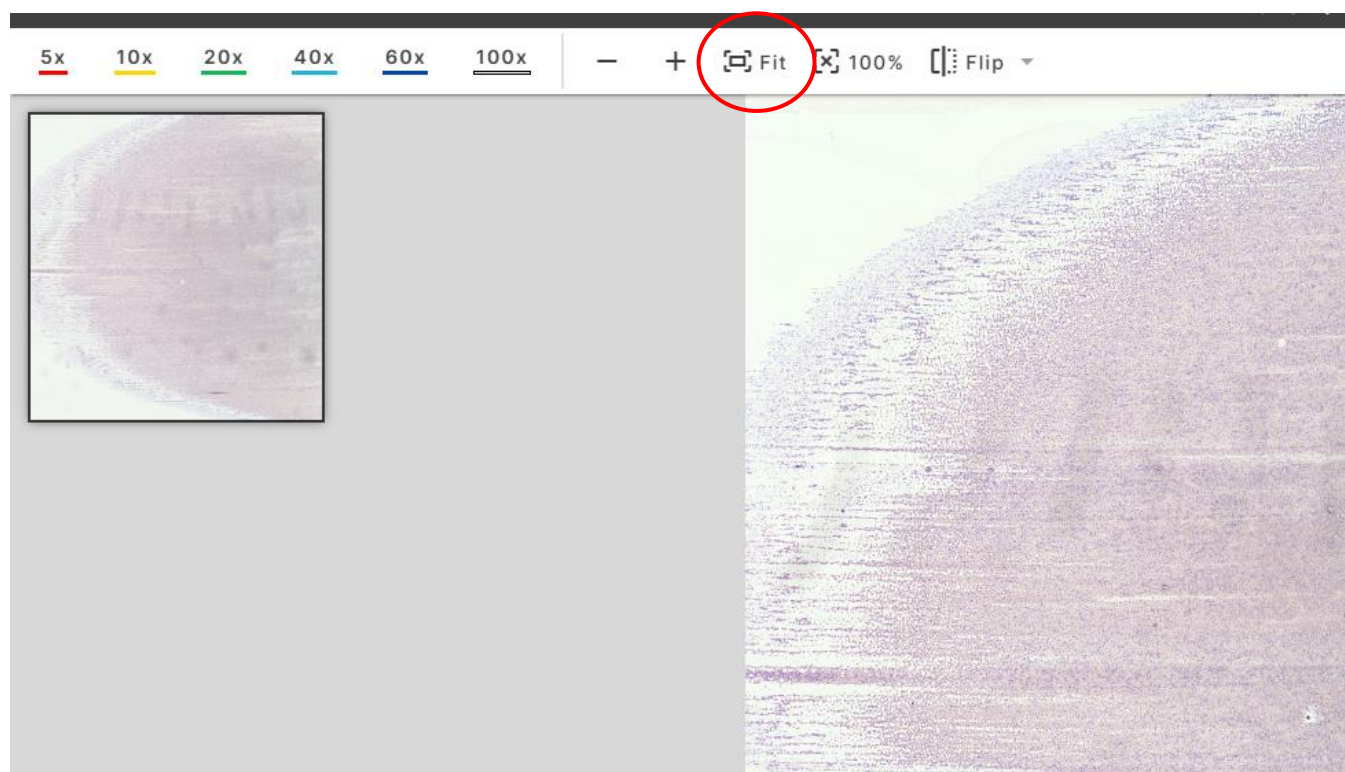

d. The 100% icon will immediately magnify the image to x100



e. The scroll wheel on the mouse will zoom in or out of the image depending on which direction is scrolled


4. The  icon will show the image as full screen. Once in full screen, clicking on the  icon will exit full screen mode and return the page to the initial view.

5. The  icon (View settings) will allow addition and removal of the minimap, the scale bar and the magnifications.

These can be toggled on and off



6. Clicking on the  icon (View Image filters) allows for adjustment of brightness, colour, contrast, sharpness etc. Please note that any changes will only affect the image on your screen and will not be saved to the image seen by other participants. Also, any adjustments can be reset using the clear all button.

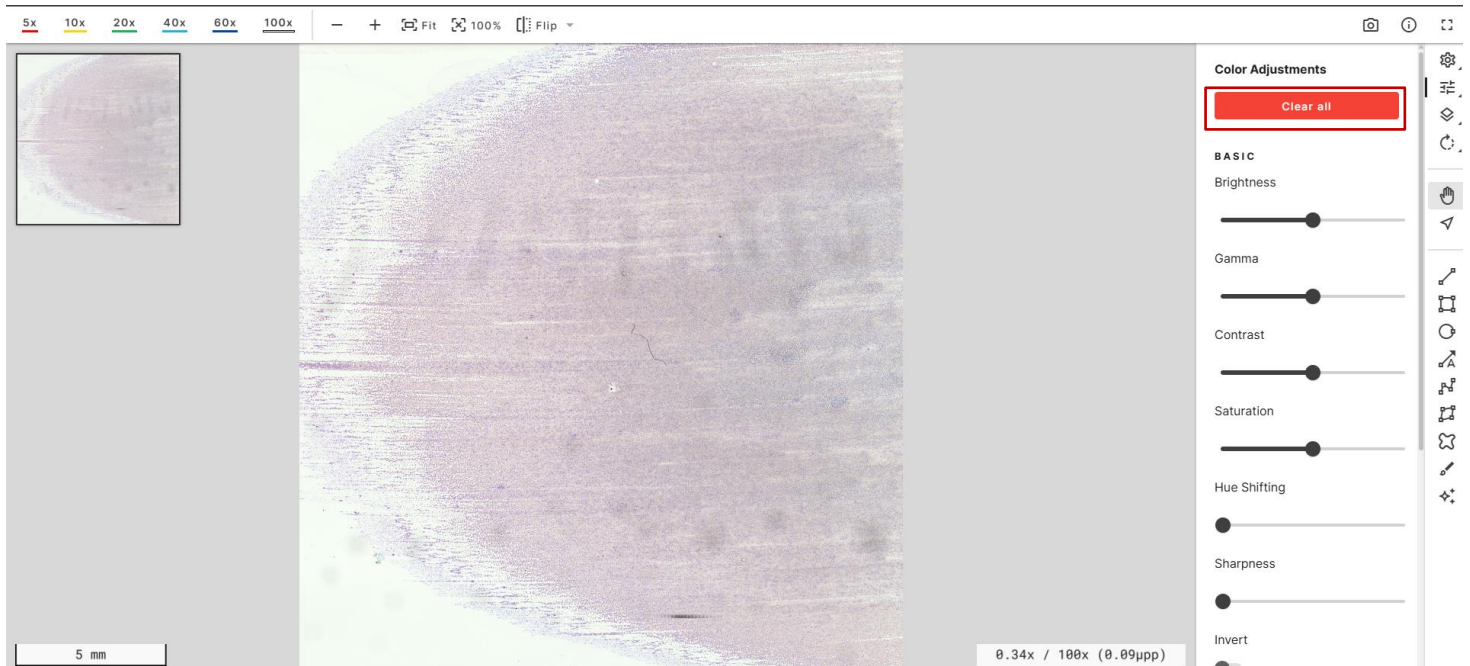

If you have any questions about accessing the image or have a problem accessing the site, please get in touch at:

UK NEQAS for Leucocyte Immunophenotyping

Telephone: (+44) 0114 267 3600

Email: admin@ukneqasli.co.uk